


KNEES

E-PRK401

- Polycentric knee
- 4-bar linkage is made of aircraft alloy (AA7075)
- Proximal attachment has A-P slide and rotation adjustment
- Suited for k2 level
- Optional adapters: Pyramid (AK) and Lotus (KD)


P/N	Pyramid Adapter	LotusAdapter
Weight Limit	125 kg / 275lbs	125 kg / 275lbs
Build Height	160 mm	160 mm
KneeCenter to Top	24.5mm	24.5mm
Product Weight	651g	672g
Material	Aluminum	Aluminum
Warranty	2 year	2 year

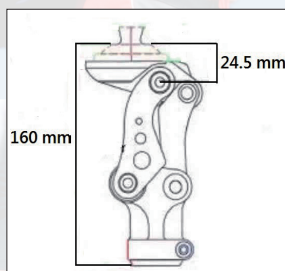

E-PRK431

- Polycentric knee joint with pneumatic swing control
- Independent flexion / extension adjustment
- All axes have built-in bearing which allow for ultra smooth walking movement
- Frame construction made of superlight aluminum alloy
- Proximal attachment has A-P slide and rotating adjustment
- Best suited for moderate k3 level
- Adjustable Extension assist
- Optional adapters: Pyramid (AK) and Lotus (KD)

P/N	Lotus Adapter	Pyramid Adapter
Weight Limit	125KG / 275lbs	125KG / 275lbs
Build Height	218mm	218mm
KneeCenter to Top	24.5mm	24.5mm
Knee Flexion	135 Degrees	135 Degrees
Product Weight	875g	897g
Material	Aluminum	Aluminum
Warranty	2 year	2 year

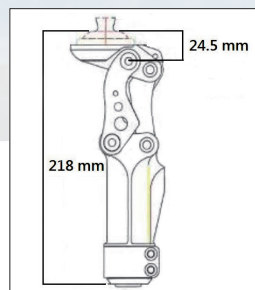

E-PRK531

- Very stable 5 Bar Knee with a large stance flexion angle (up to 12°)
- Pneumatic cylinder with independent flexion extension adjustment for swing phase
- All axes have built-in bearing which allow for ultra smooth walking movement
- Extensive alignment adjusting abilities for a smooth and natural gait pattern
- Polycentric knee with stance flexion control for maximum safety during stance
- Best suited to be used for low to high K3 level
- Practitioner adjustable from extreme stability to high stance flexion
- Adjustable extension assist
- Optional adapters: Pyramid (AK) and Lotus (KD)



	Pyramid Adapter	Lotus Adapter	
Weight Limit	125KG / 275lbs	125KG / 275lbs	
Build Height	219mm	219mm	
KneeCenter to Top	24.5mm	24.5mm	
Knee Flexion	135 Degrees	135 Degrees	
Product Weight	976g	998g	
Material	Aluminum	Aluminum	
Warranty	2 year	2 year	

LOCKS



E-PRLO-01

- Weight limit: 136 kg
- Versatile Cylindrical Body with a 4-Hole Attachment
- Serrated Pin Type Attachment
- Deep Pin Well for Easy Insertion
- Includes 2 Pins
- Variable Pin length available
- Cost Effective Pricing
- Warranty : 1 year



E-PRLO-02

- Weight limited: 114 kg
- Smooth gear mechanism
- All-in-one design
- Quiet operation
- Cost effective pricing
- Ideal for laminated sockets
- Warranty : 1 year



E-PRLO-03

- Weight limited: 100kg
- Smooth gear mechanism
- Built-in pyramid
- Screw-in modular design provides easy maintenance when required
- Cost effective pricing
- Need to use with titanium or stainless steel 3-prong or 4-prong rotating base and 4 hole rotating base
- Warranty : 1 year