

elife

enhance
life.



2020

SPORTS
MEDICINE

Europe

Austria
Belgium
Croatia
Cyprus
Denmark
France
Germany
Greece
Hungary
Italy
Macedonia
Netherlands
Poland
Russia
Serbia
Switzerland
Spain
Sweden
Turkey
UK




tel: +886-2-2923-1516
fax: +886-2-2929-0783
info@e-lifebracing.com

WORLD PARTNERS

elife is dedicated to restoring mobility and improving the quality of life for people around the world. Through our commitment to quality, innovation and product safety we deliver life enhancing product solutions that have a positive impact on a broad range of musculoskeletal conditions.

As one of the world's leading providers of orthopedic supports, e-life has built its reputation on innovative product solutions that are founded on a total commitment to quality and safety. We control every stage of the manufacturing process and apply the same exacting standards, from the sourcing of raw materials to research and development, manufacture and testing. Our commitment to quality and safety is paramount and we exceed all international quality and safety standards, including DNV, DIN, CE and FDA –

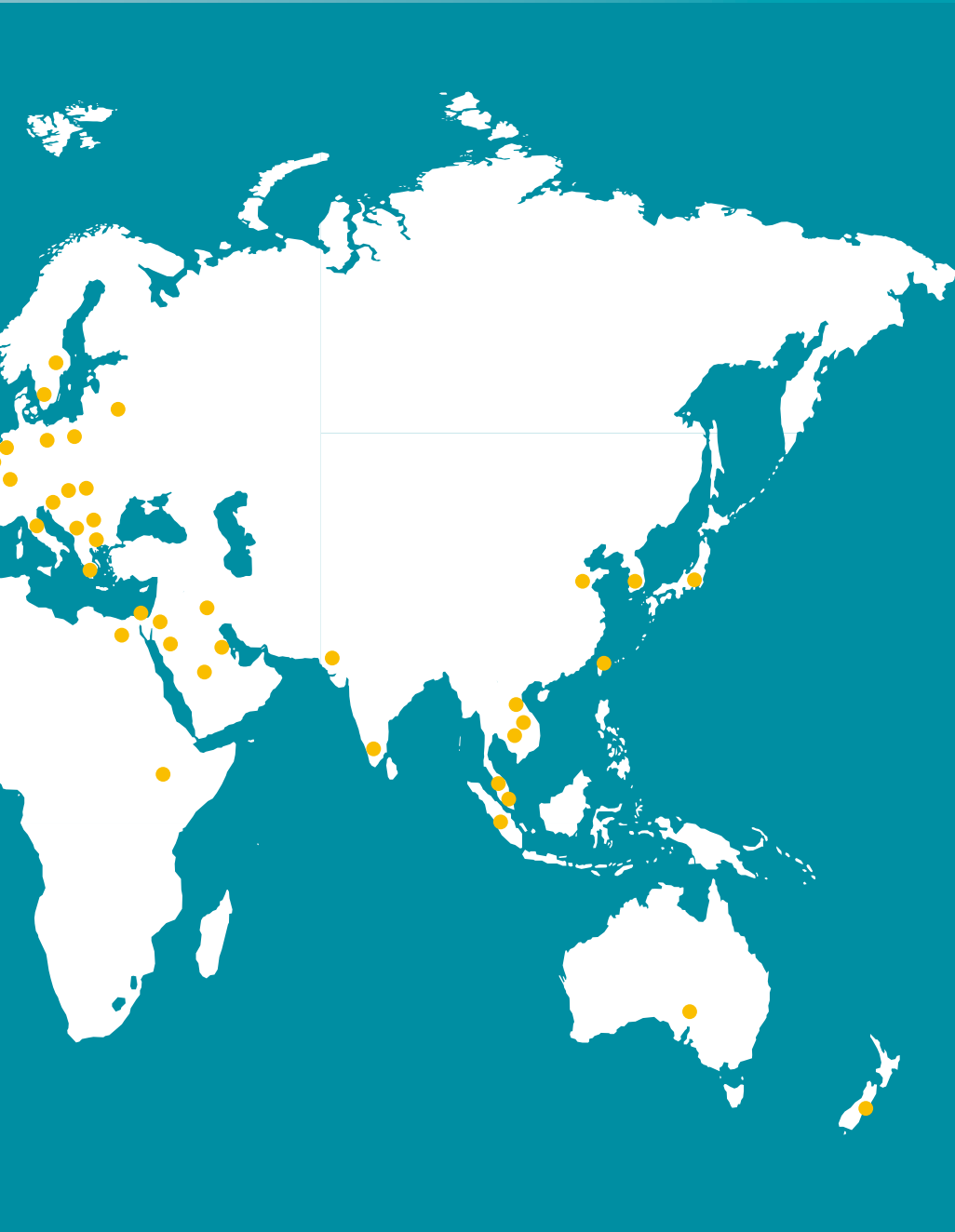

Asia Pacific

Australia
Japan
India
Indonesia
Korea
Malaysia
Singapore
Taiwan
Thailand
New Zealand
Cambodia

America

USA
Canada
Mexico
Uruguay
Colombia
Brasil

Africa

Egypt
Ethiopia

Middle East

Bahrain
Israel
Iran
Jordan
Kuwait
Lebanon
Iraq
Qatar
Saudi Arabia
UAE

we are accredited with ISO 9001 and ISO 13485. Our in-depth medical knowledge combined with our advanced product design expertise enables us to deliver functional products for optimal patient care, helping solve mobility problems and deliver a better quality of life.

Orthopedics Supportive and flexible to optimize movement and comfort. We are committed to restoring mobility and improving the quality of life for patients, with a range of braces and supports that have been developed to enhance mobility across all areas of the human body. As a world leader in medicinal orthopedics, we have the expertise and technical knowledge that ensures our products deliver maximum functionality and comfort.





CONTENT

NEW ITEMS

E-KN045 & E-KN043	8
F.A.S.T SYSTEM	11
E-WK018	13
E-WA027	14
E-HP011	16

SERIES PRODUCTS

E-KNIT SERIES	20
Q-FIT SERIES	22
D1 SERIES	26
BIO-CONTOUR SERIES	30

UNIT


KNEE	33
ANKLE	44
BACK	51
NECK	60
SHOULDER	64
ELBOW	68
WRIST	71
FOOTCARE	77
PEDIATRIC	81



Copyright © 2019 E-life International Co., Ltd. All Rights Reserved.



P.8

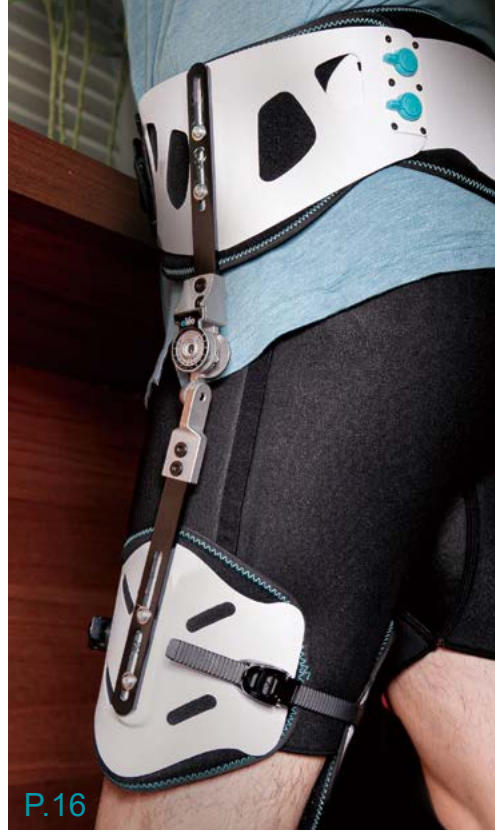

P.16



P.22



P.11





NEW ITEMS

E-life as one of the most respected orthopedic device manufacturers, we are dedicated to restoring mobility and improving the quality of life for people around the world. E-life's technology helps guarantee highest quality standards in medical brace and orthosis. Recently, many people suffer orthopedic disorders like knee pain, OA hip and spinal spondylosis, etc. Therefore, the products to correspond to above disorders are the development mainstream in the world's orthopedic device.

In the past two years, we committed to design many innovative orthoses. Continuous improvements of the technology allows for tool-free devices. It's important for us to guarantee the medical efficacy and outstanding quality of our products. We must also develop products that suit people's need. We provide new functional knee brace and hip orthosis. These products have special design for tool-free using.

E-KN045

TELESCOPIC OA KNEE

FEATURES

- The varus/valgus hinge can adjust lateral and medial shift angle to unload shear force of the knee
- Length adjustment from 16" to 17.5"
- Angular hinge allow safety control for knee motion
- Elastic outline fit thigh and shin shape and provides extremely comfort
- Lightweight, and ergonomic frame designed for superior knee support
- Unique buckle design for easy donning and doffing



Press button to adjust
strut length.

SIZE
UNIVERSAL
LEFT/RIGHT




E-KN043

COMPACT KNEE BRACE

FEATURES


- Tool-less ROM Hinge with stops cover makes sure angle stop in right place.
- ROM hinges with plug-in flexion/extension stops.
Extension limitation at:0°,10°,20°,30°,40°
Flexion limitation at:0°,10°, 20°,30°,45°,60°,75°,90°
Immobilization limitation at: 0°,10°,20°,30°
- Ultra-light weight aluminum alloy frames decrease the loading on knee joint.
- Heat molded EVA foam paddings without stitching seam and have anti-slip silicon.
- Cushioned pads with anti-slip silicon prevent slipping down.
- Labeled number straps guide the application sequence clearly.

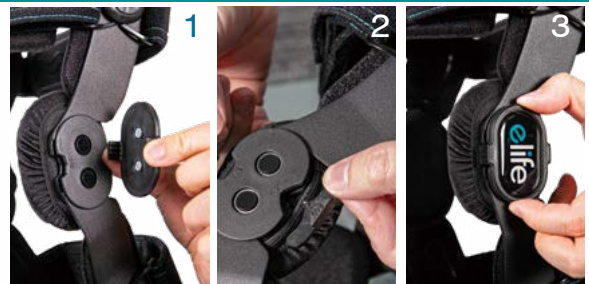

- 1.Remove magnetic cover
- 2.Plug in stops
- 3.Put the magnetic cover back

SIZE	THIGH(CM)	CALF(CM)
S	31-39	21-29
M	39-47	29-37
L	47-54	37-43
XL	54-60	43-50

LEFT /RIGHT





Inflate fast



Deflate easily


F.A.S.T SYSTEM

PATENT
PENDING

F.A.S.T series products use innovation technique to form air column which would provide support and compression. Tool-free design with embedded pump inflates air quickly and easily. Users can adjust fitting level by embedded pump. Additional inflated area design let customers feel less stress on injured tissue when inflating. F.A.S.T series products use skin friendly materials for absorbing moisture.



Air column

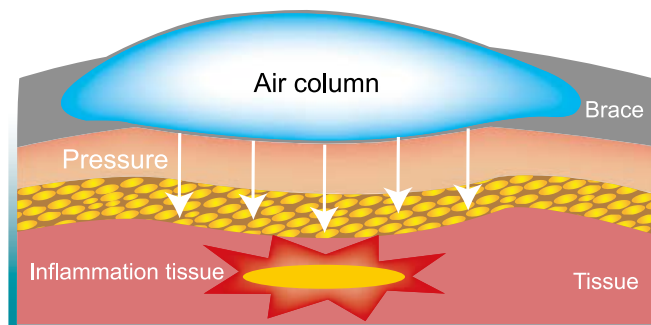

TARGETED PRESSURE TO RELIEVE PAIN



E-WR066
F.A.S.T
WRIST BRACE

FEATURES

- Adjustable straps for customized fit.
- Fit either right or left wrist.



E-WR067
F.A.S.T
WRIST SPLINT

FEATURES

- Anatomically-contoured aluminium palmar splint.
- Adjustable straps for customized fit.
- Fit either right or left wrist.



SIZE	CM	
S	13-16	
M	16-19	
L	19-22	

Suitable for E-WR066, E-WR067

E-AN066
F.A.S.T
ANKLE BRACE

FEATURES

- The injected air pressure provides extra support on the unstable ankle joint.
- Designed with elastic straps for tension adjustment.



E-KN066
F.A.S.T
KNEE BRACE

FEATURES

- The injected air pressure relieves the strain on the knee without interference with the activity
- Designed with elastic straps for tension adjustment.



SIZE	A (CM)	B (CM)	
S/M	23-28	20-25	
L/XL	28-33	25-30	
LEFT /RIGHT			

SIZE	CM	
S	33-35.5	
M	35.5-38.5	
L	38.5-43	



NEW ITEM 

E-WK018 **TRANS 2 WALKER**

FEATURES

- Double walker pumps can effortlessly and conveniently use.
- Lightweight design including full shell.
- Front panel for more complete protection and support.
- Athletic outsole design provides better shock absorption and improves anti-slip ability.

PATENT NO.
DE202019103769U1



Pull rod to let air in pump



Push rod to inflate airbag



Deflated button

SIZE	WOMEN SHOE (US)	MEN SHOE (US)
S	6-8	4-7
M	8-11	7-10
L	11-15	10-13
XL	15+	13+

E-WA027

APS COMBO BACK



FEATURES

- Flip cover design with back semi-rigid panel
- Adjustable string and roller system perfectly supports the core muscle
- Easy add accessories: front and/or side panels
- Available large or small back panel
- Universal size design fits most users

OPTIONS



Small front panel



Side panel



Large back panel



Large front panel



Additional front fabric has clearly size information. Users could extend the brace by using additional front fabric. Then, users adjust different sizes by adjustable hook and loop closure.

S	70-80 cm
M	80-90 cm
L	90-100 cm
XL	100-110 cm
2XL	110-120 cm
3XL	120-130 cm
4XL	130-140 cm

SIZE	
UNIVERSAL	

APS COMBO BACK limits motion and provides effective trunk support, and superior comfort. This brace is available in a universal size. Breathable materials and an ergonomic design, APS COMBO BACK contours comfortably to the vast majority of customers anatomies. All accessories components are easily assembled to further customize the brace to accommodate various patient conditions.



E-HP011

ADJUSTABLE HIP ORTHOSIS

PATENT
PENDING



Rotate component for adjusting waist circumference.

FEATURES

- Special component for adjusting waist circumference easily and quickly.
- Double hip hinges can set ROM limitation individually.
- One push buckle design for micro adjustment of tightness over thigh cuff.
- Rigid outer shell on waist and thigh for great stability.
- Flexion/extension ROM control from 0° to 120° in 10° stop settings, abduction - adduction adjustments are allowed.
- The universal hinge assembly can be moved to the opposite side and offers a height adjustment at the hip and thigh.
- Light-weight design to improve customers comfort and compliance.
- Breathable and washable padding liner on waist and thigh for optimal skin contact.

Double hinges could set ROM limitation individually.



One push buckle design for micro-adjustment of tightness over thigh cuff.

SIZE	Pelvic (CM)	Proximal thigh (CM)	
S/M	66-96	40-65	
L/XL	96-127		

ADJUSTABLE HIP ORTHOSIS offers quality and hip motion restriction for pre-, post- and non-operative patients. It uses pended patent component to adjust hip waist circumference. Light-weight design to improve customers comfort and compliance without compromising support.



SERIES

As we develop more kind of products in recent years, we are rich of complete product line. Therefore, there are four series type of products in E-life. They are classified by material characteristics, similar mechanism or their intended use. We have E-knit series, Q-Fit system series, D1 series and Bio-Contour braces series.

E-knit series products are ideal for symptomatic support while engaging in active mobility. Q-Fit system series products use BOA® system to don and doff. D1series products mold from two different materials. Bio-Contour braces series products have innovative 3D flat knitting technology. Bio-Contour braces series products are inserted silicon pad or supportive stay.

E-KNIT SERIES


E-knit series products are made of circular-knitted, four-way stretch elastic material for comfortable fit. E-life knitting technique produces three-dimensional pattern which is located around the peripheral joints or major muscle groups. 3D patterns provide constant and stable stimulation that improve motion control by stimulating proprioception receptors. E-knit series products are ideal for symptomatic support while engaging in active mobility.



E-KN202
GENU-KNIT
BRACE




SIZE	CM
S	28-34
M	34-41
L	41-48
XL	48-55
XXL	55-63



E-AN202
MALLEO-KNIT
BRACE




SIZE	CM
S	15-20
M	20-25
L	25-30
XL	30-35

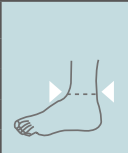


E-TH202
FEMORIS-KNIT
BRACE




SIZE	CM
S	44-49
M	49-54
L	54-59
XL	59-65



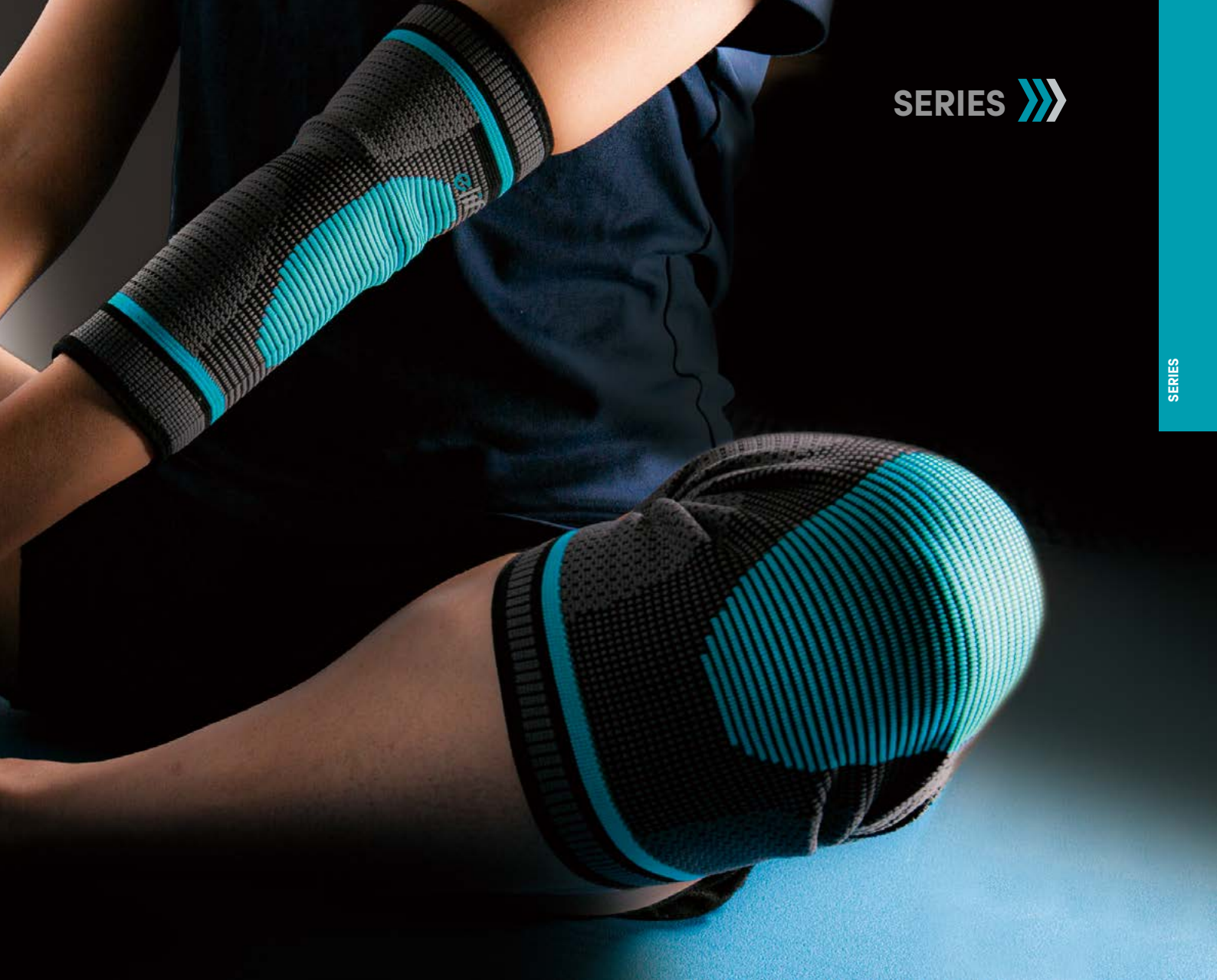


E-SN202
CALF-KNIT
BRACE



SIZE	CM	
S	30-34	
M	34-38	
L	38-42	
XL	42-46	


E-EL202
EPI-KNIT BRACE



SIZE	CM	
S	23-26	
M	26-31	
L	31-37	
XL	37-43	

E-WR204
MANU-KNIT
BRACE

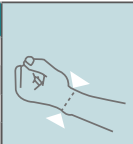

SIZE	CM	
S	13-15	
M	15-18	
L	18-20	
XL	20-22	

BOA®


Q-FIT SERIES

The Q-Fit system is developed by conception of easy and quick donning and doffing. More than that it brings users more benefits. Q-Fit system replaced the traditional lace function so that users can be convenient to use and save more time and energy during taking on or off the brace.

During activities, users are able to enhance effectiveness of performance and excellent compliance.

E-WAR021A PREMIUM Q-FIT BACK ARMOUR

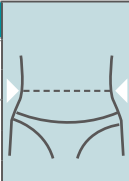


SIZE	CM	BACK PANEL	
S	75-85	S	
M	85-95		
L	95-105	L	
XL	105-115		
XXL	115-125		

E-WAR110 Q-FIT LOW PROFILE BACK BRACE

FEATURES

- Low profile design with 4 stays used to support L1-L5
- Easy to adjust intra-abdominal pressure when changing position from sitting.
- Pocket handle for convenient applying

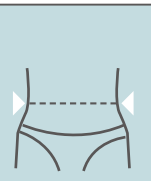

SIZE	CM	
S	56-71	
M	69-84	
L	81-96	
XL	94-109	
XXL	107-122	



E-WAR021 Q-FIT BACK ARMOUR

FEATURES

- The low profile posterior panel, anterior panel and side panels all are connected via a Velcro closure, and can be removed and reattached easily. (E-WAR021A)
- The Q-Fit system allows you to do quick, one-handed control your comfort for standing or sitting.
- Thin, lightweight, breathable mesh material combines with contouring panels for increased support.




E-HP010 FAN-LITE HIP ORTHOSIS

FEATURES

- Q-fit system with single dial front closure
- Semi-universal-sizing-circumference tape allows for accurate patient sizing.
- Small joint head and lightweight design with 0-120 flexion and extension adjustment
- Tool-free length adjustment design allows to position thigh pad and thigh cuff

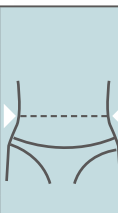

	SIZE	CM	
Regular	S	75-85	
	M	85-95	
	L	95-105	
Large	XL	105-115	
	2XL	115-125	
	3XL	125-135	



E-WRR052 Q-FIT WRIST SPLINT

FEATURES

- Anatomically-contoured aluminum palmar splint
- Immobilization of the wrist
- Comfortable, breathable, skin-friendly material

SIZE	CM	
S	13-15	
M	15-18	
L	18-20	
XL	20-22	
LEFT /RIGHT		

E-WRR058 Q-FIT THUMB AND WRIST SPLINT

FEATURES

- Anatomically-contoured aluminum palmar splint
- Immobilization of the wrist and thumb
- V-notch at fifth MP for palmar spread
- Breathable open cell foam material

SIZE	CM	
S	13-15	
M	15-18	
L	18-20	
XL	20-22	
LEFT /RIGHT		



E-ELR006 Q-FIT ELBOW BRACE

FEATURES

- The SBR pad provides pain relief at the muscle insertion point
- Extra elastic band fixes brace in position for easy donning and doffing

SIZE	CM	
S	23-25.5	
M	25.5-28	
L	28-30.5	
XL	30.5-33	




E-ANR043 Q-FIT ANKLE BRACE

FEATURES

- Bilateral plastic stays for additional support
- Extra non-elastic cross-stirrup straps provide maximal stability
- Made of duration and tough nylon fabric



SIZE	CM
S	15-20
M	20-25
L	25-30
XL	30-35




E-ANR045 Q-FIT ANKLE STABILIZER

FEATURES

- Y-Strap stabilizes the ankle and prevent it from rolling inwards by Q-Fit system adjusting tension on the straps.
- Bilateral anatomically PE plates provide great support
- Easy rear foot entry allows for quick on-off



SIZE	CM
S	15-20
M	20-25
L	25-30
LEFT /RIGHT	



E-AFK002 Q-FIT CARBON AFO

FEATURES

- Carbon fiber strut with anatomical shape for superior fit
- Easy to don and doff with Q-Fit system
- Trimmed footplate for fitting most shoes
- Designed specifically for right or left




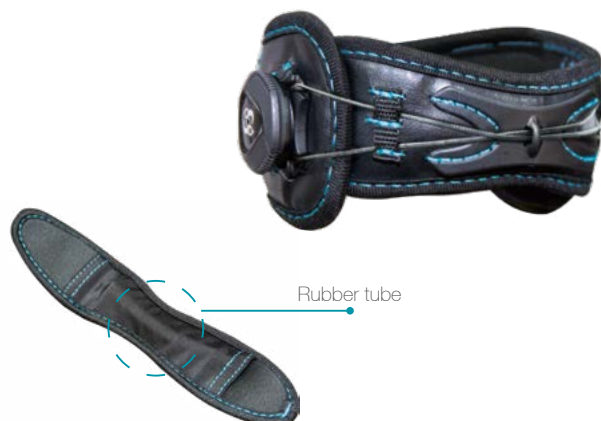
SIZE	SHOE SIZE (EU)	HEIGHT(CM)
S	35-37	32
M	37-39	34
L	39-41	36
XL	41-44	38
LEFT /RIGHT		

E-KNR005 Q-FIT PATELLAR STRAP

FEATURES

- Click buckle designed is more convenient for using
- The rubber tube provides the patellar tendon support to release pain
- Breathable material for great comfort

SIZE	CM
S	30-34
M	34-38
L	38-42
XL	42-46
XXL	46-50

D1 SERIES

D1 series is an innovative manufacturing process used to produce complicated molded parts from two different materials. D1 series molding optimizes to create powerful combination of hard and soft materials and enhanced tactile and cosmetic product features.

E-WK016

D1 WALKER

FEATURES

- Double-injection design with full shell and fancy style provides comfort for calf area
- Double pumps provide individual compression to decrease the swelling
- Shock absorbing sole reduces heel-strike impact force
- Toe guard for more complete protection and support
- Athletic outsole design provides better shock absorption and improves anti-slip ability
- 11" and 17" available



Double pumps

WK016-11"



E-WK017

D1 WALKER WITH AIR CUSHIONED OUTSOLE

- Air-cushioned outsole design provides excellent shock absorption

SIZE	WOMEN SHOE	MEN SHOE
S	6-8	4-7
M	8-11	7-10
L	11-15	10-13
XL	15+	13+

Suitable for E-WK016, E-WK017



Air-cushion



These options are available for elife walkers.

FEATURES

- The varus/valgus hinge can adjust lateral and medial shift angle to unload shear force of the knee (E-KN041)
- Angular hinge allow safety control for knee motion (E-KN041)
- Medially applied frame corrects for either medical or lateral compartment
- Lightweight, and ergonomic frame designed for superior knee support
- Unique buckle design for easy donning and doffing
- Breathable liner provides optimal texture

SIZE
UNIVERSAL
LEFT/RIGHT

**E-KN041
OA KNEE SHELTER**

**E-KN040
OA KNEE
STANDARD**




E-WR009 CMC GUIDER

FEATURES

- Ergonomic frame provides superior support and excellent comfort of CMC joint
- Inner metal stay offers extra support and stability of thumb
- Easy to apply and adjust with simple straps

SIZE	CM
S	16-19
M	19-23
L	23-26
LEFT /RIGHT	




E-AN065 D1 STIRRUP BRACE

FEATURES

- Anatomically-contoured plastic shell
- Inflatable air bladders allow user to adjust compression pressure to reduce swelling
- Simple to apply and adjust by two circumferential straps
- Easy adjustable heel strap

SIZE
UNIVERSAL



E-AN064 DYNAMIC ANKLE SHELTER

FEATURES

- Adjustable module is designed for customize fit
- Hinges allow free plantar motion and prevents damaged ankle torsion
- Bilateral soft pad eliminates joint pressure
- Moldable foot plate enhances ankle control and comfort

SIZE	WOMEN SHOE	MEN SHOE
S/M	8-11	7-10
L/XL	11-16	10-15



BIO-CONTOUR SERIES

Innovative 3D flat knitting technology makes braces provide optimal body fit. Bio-contour series perfectly provides support over the joints, muscles and soft tissue. Bio-contour braces not only enhance the injury joints and soft tissue recovery ability but also prevent the injure happened during activities.




E-WA401 BIO-CONTOUR BACK BRACE WITH STRAP

FEATURES

- 3 plastic stays provide extra low back support
- Bilateral compression bandages provide intra-abdominal pressure
- Hand pocket design for easy donning and doffing.

SIZE	CM
S	68-78
M	78-88
L	88-98
XL	98-108
XXL	108-118

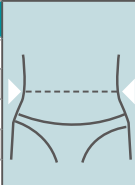



E-WA450 BIO-CONTOUR BACK BRACE

FEATURES

- Triangular silicone pad provides great low back support.
- 3 plastic stays offer extra back stability.
- Hand pocket design for easy donning and doffing.


SIZE	CM
S	68-78
M	78-88
L	88-98
XL	98-108
XXL	108-118

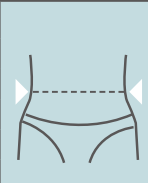


E-KN450
BIO-CONTOUR KNEE BRACE



SIZE	CM
S	33-36
M	36-38
L	38-41
XL	41-46



E-KN451
BIO-CONTOUR HINGED KNEE BRACE

FEATURES

- 3 dimension knitting creates an anatomical fit.
- Collateral hinges prevents valgus or varus deviation and for hypertension control.(E-KN451)
- Silicone pad maintains patella in right position.
- Straps for individual compression adjustment. (E-KN451)




E-AN450
BIO-CONTOUR ANKLE BRACE

FEATURES

- 3 dimension knitting creates an anatomical fit
- Silicone pad maintains angle joint in right position
- Designed with elastic straps for tension adjustment. (E-AN451)



SIZE	CM
S	15-20
M	20-25
L	25-30
XL	30-35



E-AN451
BIO-CONTOUR ANKLE BRACE WITH STRAP




E-SH450 BIO-CONTOUR SHOULDER BRACE

FEATURES

- Special silicone pad improves muscle activation
- Dual adjustable straps for easy using and fit both sides
- Axillary soft pad design for comfortable wearing



SIZE	CM
S	25-28
M	28-31
L	31-34
XL	34-37




E-EL450 BIO-CONTOUR ELBOW BRACE

FEATURES

- Two silicone pads provide intermittent massage to targeted areas.
- The 3D knitted material of the elbow brace secures it in place



SIZE	CM
S	21-23
M	23-25
L	25-28
XL	28-31




E-WR450 BIO-CONTOUR WRIST BRACE

FEATURES

- Provide great support and compression.
- Silicone pad prevents over-stretch on ulnar side




SIZE	CM
S	13-15
M	15-18
L	18-20
XL	20-22

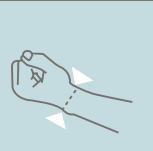

KNEE

POST OP KNEE BRACE

Post-Operative Knee Brace is designed to provide protected, controlled range-of-motion for patients recovering from knee surgery or those who have knee injuries or instabilities as following indications:

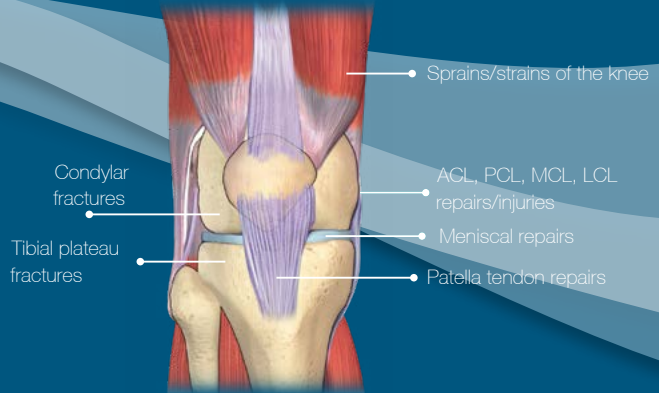

E-KN099 DELUXE KNEE BRACE

FEATURES

- Streamlined durable aluminum hinge provides maximum protection but with minimum weight
- Simply press and rotate hinge with lock
 - Extension limitation from 0° ~ 105°
 - Flexion limitation from 0° ~ 120°
 - Lock from 0° ~ 90°
 - Adjustable in every 15 degrees
- "Pull pin design" enables easy release for slide-to-size length adjustments.



MALLEABLE ALUMINUM CUFFS



SIMPLY PRESS HINGE LOCK



ANTI-SLIP PAD



INFO	LENGTH (CM)	CIRCUM. (CM)	SIZE
E-KN099	APPROX 48-61		UNIVERSAL 
E-KN099B	LONG 57 / SHORT 44	CALF MAX.48 THIGH MAX. 61	
E-KN099C	APPROX 48-61		

E-KN099B
FULL FOAM KNEE BRACE



E-KN099C
COOL FOAM KNEE BRACE



CUTTABLE STRUTS



LENGTH ADJUSTMENT



COOL FOAM



E-KN097B

POST OP KNEE BRACE

FEATURES

- Excellent collateral stabilization via long side-bars prevents valgus or varus deviation.
- Adjustable ROM by simple setting hinges
- Telescoping feature for tool-free length adjustment
- Simply rotate and press-lock settings
 - Extension limitation from 0° ~ 60°
 - Flexion limitation from 0° ~ 120°
 - Immobilization limitation from 0° ~ 60°
- Adjustable in every 15 degrees



E-KN098

SETTER KNEE BRACE

FEATURES

- "Unique lock design" enables easy release for slide-to-size length adjustments.
- Simply press and rotate hinge with lock
 - Extension limitation from 0° ~ 60°
 - Flexion limitation from 0° ~ 120°
 - Lock from 0° ~ 45
 - Adjustable in every 15 degrees
- Easily apply on and take off by side buckles
- Length adjustable from 19" to 24"


SIZE	LENGTH (CM)	CIRCUM. (CM)	
UNIVERSAL	APPROX 48-61	CALF MAX.48 THIGH MAX. 61	

Suitable for E-KN097B, E-KN098

OA KNEE

In people who have osteoarthritis (OA) knee, the cartilage wears away, leading to an overgrowth of bone underneath. The cartilage becomes rough and then breaks down.

When the cartilage breaks down, it causes pain, swelling, and problems with moving the joint. These symptoms can worsen over time.

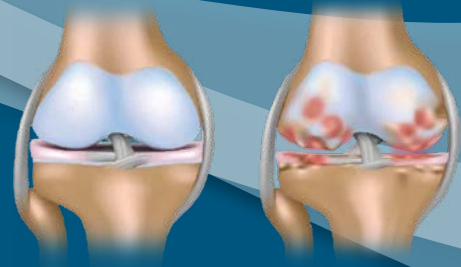

Normal knee

Osteoarthritis (OA) knee


E-KN042

DAILY OA KNEE BRACE

FEATURES


- Provide offloading of the medial or lateral compartment for mild to moderate osteoarthritis.
- The single-sided low profile hinge encourages patient compliance
- “U” Buttress controls the patella in a more natural position.
- Open anterior thigh section for easy donning and doffing
- Comfortable, breathable material provides focused compression around the knee joint to address swelling and promote blood flow

OFFLOADING MECHANISM

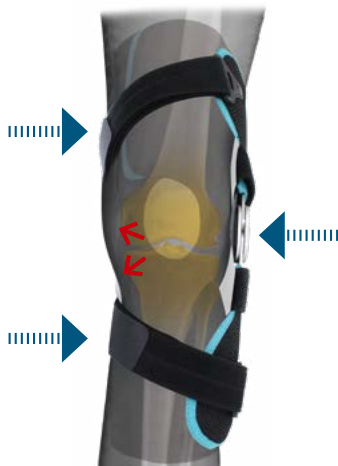

3 POINTS OF LEVERAGE

SIZE	THIGH(CM)	CALF(CM)
S	30-37	25-30
M	37-44	30-35
L	44-52	35-42
XL	52-60	42-50

LEFT /RIGHT 





E-KN036

**UNLIMIT S1 OA
KNEE BRACE**

FEATURES

- Ultra-light weight, single upright with valgus and varus adjustment to unload
- Laterally applied frame corrects for either medical or lateral compartment
- 3 points of leverage for offloading.
- Buckle design eliminates Velcro adjustment when every donning and doffing.
- One size fits most users and allows ergonomically fitting
- ROM hinges with plug-in flex/ext stops

SIZE
UNIVERSAL
LEFT/RIGHT

E-KN038


**CONQUER OA
KNEE BRACE**

FEATURES

- Dual uprights with valgus and varus adjustment for offloading
- Double upright provides functional support for those patients with ligament instability
- Support the knee joint and share the loading of affected joint
- ROM hinges with plug-in flex/ext stops
- Full shell design provide well stability
- Could be used at home or for gentle sports




SIZE	THIGH(CM)	CALF(CM)
S	39-47	29-37
M	47-54	37-43
L	54-60	43-50
XL	60-67	50-57
LEFT /RIGHT		



FUNCTIONAL KNEE BRACE

Functional Knee Brace can effectively provide great stability for tears, injuries, or reconstruction of knee ligaments as followings:

- Anterior cruciate ligament (ACL)
- Posterior cruciate ligament (PCL)
- Lateral collateral ligament (LCL)
- Medial collateral ligament (MCL)

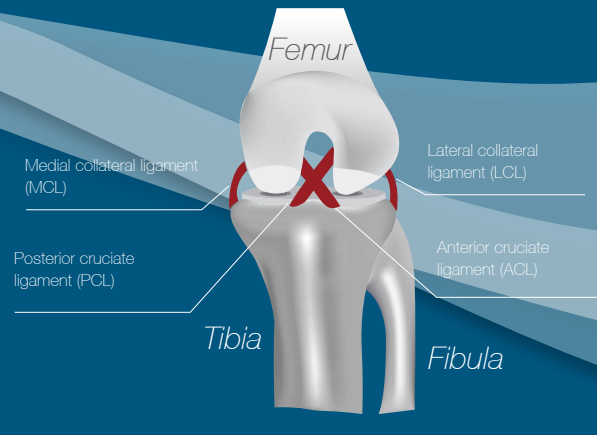

E-KN034 FUNCTIONAL KNEE BRACE

FEATURES

- 4-leverage design for optimal stabilization
- Prevent hyperextension via adjustable hinge ROM
- ROM hinges with plug-in flex/ext stops
- Extension limitation at: 0°, 10°, 20°, 30°, 40°
Flexion limitation at: 0°, 10°, 20°, 30°, 45°, 60°, 75°, 90°
Immobilization limitation at: 0°, 10°, 20°, 30°
- Silicone-grip strap pads for anti-migration
- Concentric condyle pads correspond to bilateral femoral condyles.

1. Point of leverage
Femur is anchored by a cuff placed at the top of the anterior thigh.




2. Point of leverage
A strap across the back of the lower thigh pushes the femur anteriorly

3. Point of leverage
A strap placed on the anterior tibia applies an active constant load to prevent anterior tibial translation

4. Point of leverage
Tibia is anchored by a cuff placed bottom of the posterior calf

SIZE	THIGH(CM)	CALF(CM)
S	31-39	21-29
M	39-47	29-37
L	47-54	37-43
XL	54-60	43-50
LEFT /RIGHT		




E-KN002 KNEE BRACE WITH PATELLAR OPEN

FEATURES

- 1/8" thick neoprene with nylon lining
- Increase blood circulation through warm effect
- Open patellar design
- Effectively stabilizes the kneecap to eliminate anterior knee pain

SIZE	CM
S	30-34
M	34-38
L	38-42
XL	42-46
XXL	46-50

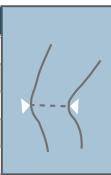



E-KN004 OPEN PATELLAR KNEE BRACE

FEATURES

- 1/8" thick neoprene with nylon lining
- Open patellar design
- Popliteal cut-out design for comfort fit
- Wrap around design provides an easy on and off with a custom fit.

SIZE
UNIVERSAL





E-KN005 PATELLAR STRAP

FEATURES

- Ideal for patellar tendinitis, jumper knee
- Reinforce and protect the patellar tendon
- Simple fitting & time-saving application via hook and loop closure

SIZE
UNIVERSAL

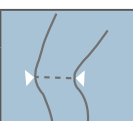





E-KN006 PATELLAR TENDON BRACE

FEATURES

- 1/8" thick neoprene with nylon lining
- Open patellar design
- Effectively stabilize the kneecap to eliminate anterior knee pain
- Wrap around design provides an easy on and off with a custom fit.

SIZE	CM
S	30-34
M	34-38
L	38-42
XL	42-46
XXL	46-50

Suitable for
E-KN006, E-KN051, E-KN052



E-KN051 HINGED KNEE SUPPORT

FEATURES

- 3/16" thick neoprene with nylon lining
- Open padded buttress supports the patella
- Lightweight hinges hidden in pockets offer superior stability
- Upper and lower fasteners offer further fit



E-KN052 HINGED KNEE SUPPORT

FEATURES

- 3/16" thick neoprene with nylon lining
- Open padded buttress supports the patella
- Lightweight hinges hidden in pockets offer superior stability
- Upper and lower fasteners offer further fit



E-KN056 UNI FIT HINGED KNEE SUPPORT

FEATURES

- Hinges hidden in pockets offer superior stability.
- Wrap around design provides an easy on and off with a custom fit.
- Upper and lower fasteners offer further fit.



E-KN057 UNI FIT ROM KNEE SUPPORT

FEATURES

- Low-profile polycentric hinges with adjustable ROM limitation.
 - Extension limitation at: 0°, 15°, 30°, 60°, 90°
 - Flexion limitation at: 0°, 15°, 30°, 60°, 90°
 - Immobilization limitation at: 0°, 15°, 30°, 60°, 90°
- Wrap around design provides an easy on and off with a custom fit.
- Upper and lower fasteners offer further fit.



SIZE	
UNIVERSAL	

Suitable for :E-KN056, E-KN057

E-KN058 WRAP AROUND HINGED KNEE SUPPORT

FEATURES

- 3/16" thick neoprene with nylon lining.
- Hinges hidden in pockets offer superior stability.
- Wrap around design provides an easy on and off with a custom fit.
- Upper and lower fasteners offer further fit.



SIZE	CM	
S	30-34	
M	34-38	
L	38-42	
XL	42-46	
XXL	46-50	


E-KN070

STABILIZED KNEE SUPPORT

FEATURES

- 1/8" thick neoprene with nylon lining
- Bilateral spiral steel stays prevent valgus or varus deviation.
- Open padded buttress supports the patella
- Upper and lower coronal fasteners offer further fit

SIZE	CM
S	30-34
M	34-38
L	38-42
XL	42-46
XXL	46-50

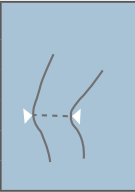


E-KN081

**KNEE
IMMOBILIZER**

FEATURES

- Three piece design for precise fit
- Removable medial, lateral and posterior stays
- Comfortable, washable, skin-friendly terry toweling cloth material
- Velcro straps with loop closures for simple fitting
- Additional patellar strapping(E-KN081A)
- Ideal for post-operative immobilization, conservative treatment of ligament injury




E-KN081A

**KNEE
IMMOBILIZER
WITH PATELLAR
STRAP**



SIZE	CM	LENGTH
S	30-34	12", 14", 16", 18", 20", 22", 24" available
M	34-38	
L	38-42	
XL	42-46	
XXL	46-50	

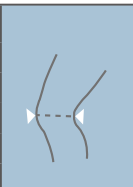



ANKLE

FEATURES

- Lightweight and low profile aluminum side struts
- Rocker bottom design assists gait
- Decreased heel-strike force on calcaneus
- Simple to apply and adjust by circumferential straps
- Dial lock hinges with adjustable ROM limitation (E-WK002)
 - Dorsiflexion limitation at: 0°, 7.5°, 15°, 22.5°, 30°, 37.5°, 45°
 - Plantarflexion limitation at: 0°, 7.5°, 15°, 22.5°, 30°, 37.5°, 45°
 - Immobilization limitation at: 0°, 7.5°, 15°, 22.5°, 30°, 37.5°, 45°
- Inflatable aircells are available for further support while reducing swelling (E-WK005)
- Available in 11" & 17" height (E-WK001/ E-WK005)
- Available in 14" & 17" height (E-WK002)

**E-WK001
FIXED WALKER**



SIZE	WOMEN SHOE (US)	MEN SHOE (US)
XS	4.5	UP TO 4
S	5-8.5	4.5-7
M	9-11	7.5-10
L	11.5-13.5	10.5-12.5
XL	13.5+	12.5+

Suitable for E-WK001, E-WK002, E-WK005

**E-WK002
ROM WALKER**



**E-WK005
PNEUMATIC FIXED WALKER**



FEATURES

- Lightweight and contoured strut design provides better support and fit
- Front and back panels for superior protection
- Rocker bottom design assists gait
- Decreased heel-strike force on calcaneus
- Simple to apply and adjust by circumferential straps
- Available in 11" & 17" height
- Inflatable aircells are available for further support (WK009/ WK009D)
- Diabetic Walker System includes: (WK009D)
 - Diabetic Walker
 - Insole WK009D
 - Stocking

SIZE	WOMEN SHOE (US)	MEN SHOE (US)
S	6-8	4.5-7
M	8.5-11.5	7.5-10
L	11.5-13.5	10.5-12.5
XL	13.5+	12.5+

Suitable for E-WK008, E-WK009, E-WK009D

E-WK008
PRO SHIELD WALKER



E-WK009
PRO SHIELD PNEUMATIC WALKER



E-WK009D
PRO SHIELD DIABETIC WALKER



E-WK014

JEEP WALKER

FEATURES

- Rubber outsole provide superior grip and abrasion resistant.
- Low profile footbed height design that practically removes leg length discrepancy.
- Ultra-light weight with anatomical frame.
- Rocker bottom encourages natural gait.
- Available in 11" & 17" height.



WK014-11"



E-WK015

PNEUMATIC JEEP WALKER

FEATURES

- Rubber outsole provide superior grip and abrasion resistant.
- Low profile footbed height design that practically removes leg length discrepancy.
- Inflatable aircells are available for further support while reducing swelling.
- Ultra-light weight with anatomical frame.
- Rocker bottom encourages natural gait.
- Available in 11" & 17" height.



WK015-11"



SIZE	WOMEN SHOE (US)	MEN SHOE (US)
XS	<5	<4.5
S	5-8.5	4.5-7
M	9-11	7.5-10
L	11.5-13.5	10.5-12.5
XL	13.5+	12.5+

Suitable for E-WK014, E-WK015

E-NS001

PLANTAR FASCIITIS STRETCH SPLINT



FEATURES

- Prefabricated splint
- Adjustable dorsiflexion straps
- Soft velour-lined foam padding for optimal comfort
- Minimal house-hold ambulation is possible via pawl surface under the sole of orthosis
- Indications of use: plantar fasciitis, achilles tendonitis, contractures and heel spurs

SIZE	WOMEN SHOE (US)	MEN SHOE (US)
S	Up to 7.5	Up to 7
M	8-10.5	7.5-10
L	11+	10.5+

E-NS003

DORSAL NIGHT SPLINT



FEATURES

- Speed recovery from plantar fasciitis, achilles tendonitis and overuse injuries
- Comfortably stretches the calf and soft tissues of the foot
- Lightweight and breathable soft padding material for comfort
- Open heel to hold heel in place and prevent migration during the nighttime

SIZE	WOMEN SHOE (US)	MEN SHOE (US)
S / M	6-10	5-9
L / XL	10.5+	9.5-14

E-AFK001

X1 CARBON AFO

FEATURES

- Carbon fiber design is light to provide excellent dynamic motion
- Easy donning and doffing with Velcro strap
- Trimmed foot part for fitting most shoes
- Anatomical shaped posterior strut for excellent fit

SIZE	SHOE SIZE (EU)	HEIGHT(CM)
S	35-37	32
M	37-39	34
L	39-41	36
XL	41-44	38
LEFT /RIGHT		



E-AF002

AFO LEAF SPRING

FEATURES

- Made of durable PP material
- Improve toe clearance in normal gait
- Designed specifically for right or left

SIZE	SHOE SIZE (EU)	HEIGHT(CM)
S	35-37	30.5
M	37-39	32
L	39-41	34
XL	41-44	39
LEFT /RIGHT		




E-AN043B

ANKLE BRACE WITH STRAP

FEATURES

- Made of durable and tough nylon fabric
- Extra non-elastic cross-stirrup straps provide maximal stability
- Bilateral plastic stays for additional support
- Velcro lace closure design for convenient application

SIZE	CM
S	15-20
M	20-25
L	25-30
XL	30-35





E-AN051

LACED ANKLE BRACE

FEATURES

- Made of durable canvas material
- Anatomically-contoured wide aluminum ankle stays
- Lace-up design

SIZE	CM
S	15-20
M	20-25
L	25-30
XL	30-35

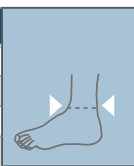


E-AN055B

AIR ANKLE STIRRUP BRACE

FEATURES

- Inflatable air bladders allow user to adjust compression pressure to reduce swelling
- Rigid contoured outer shells control inversion or eversion
- Simple to apply and adjust by two circumferential straps

SIZE
UNIVERSAL




E-AN060

AIR / GEL ANKLE BRACE

FEATURES

- Combined air and gel bladders for cold therapy and even-distributed pressure
- Simple to apply and adjust by two circumferential straps
- Easy adjustable heel strap

SIZE
UNIVERSAL





BACK 

BACK

BACK

ROLLER SPEED SYSTEM


The system is created by conception of saving energy and convenient adjust abdominal pressure. Two sides pulling provides controlled compression and also prevents brace to migration.

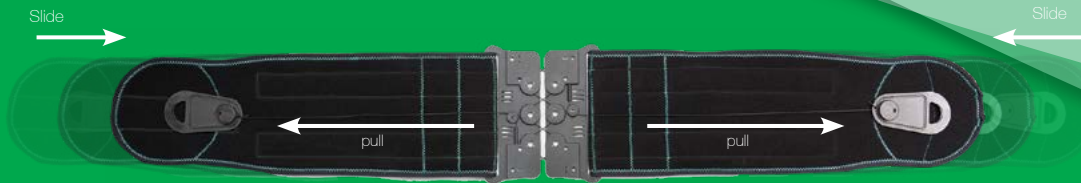

UNIVERSAL SIZE DESIGN

5 sizes in one, users can adjust sizes easily by hook and loop. Suggestion waist circumference:

S	70-85 cm
M	85-95 cm
L	95-105 cm
XL	105-115 cm
XXL	115-125 cm

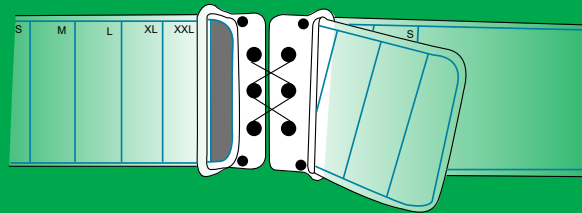

E-WA021 LSO PANEL

FEATURES

- Wrap around design with back, abdominal panel and front panel from stomach to clavicle (WA023)
- Sternal Pad Kit available for maximum flexion control (WA023)
- Anatomically-shaped rigid panel fits thoracic and lumbar curvature (WA023)



E-WA023 FRONT & TLISO PANEL






E-WA003 BACK SUPPORT

FEATURES

- Plastic splints provide extra support
- Effective vented tension side-pulls provide superior support
- Made of high quality vented elastic strapping material

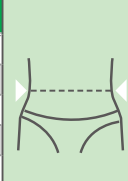
SIZE	CM	
S	56-71	
M	69-84	
L	81-96	
XL	94-109	
XXL	107-122	



E-WA004 INDUSTRIAL BACK SUPPORT

FEATURES

- Effective tension side-pulls provide superior support
- Anti-migrated textured webbing and suspension system
- Made of highly vented mesh material

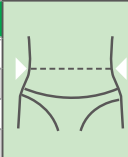SIZE	CM	
S	70-81	
M	81-96	
L	96-119	
XL	119-142	

E-WA011 FIT LUMBAR SUPPORT

FEATURES

- 4 semi-rigid steel stays support back muscles
- Made of high quality vented elastic strapping material.
- Finger pockets design with strapping system can be fitted easily without effort.




E-WA012 LUMBAFIT BACK BRACE

FEATURES

- 4 semi-rigid steel stays support back muscles
- Made of high density elastic fabric, this lightweight waist brace is soft and breathable
- With double elastic pressure straps
- Finger pockets design with strapping system for easy fitting

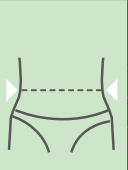

SIZE	CM	
S	56-71	
M	69-84	
L	81-96	
XL	94-109	
XXL	107-122	

Suitable for E-WA011, E-WA012

E-WA024


VENTILATED LUMBARFIT BRACE

FEATURES

- Anatomically-contoured arc stays fit lumbar curve and support lower back muscles.
- Ventilated holes provide a comfortable fit and reduce heat and moisture.
- Extra side straps for fixing V-shaped belts add convenience of wearing.
- Buttock curve contour and anatomical shape design.



SIZE	CM
S	65-80
M	75-90
L	85-100
XL	95-110
XXL	105-120

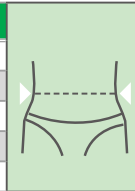

E-WA025


LIGHT LUMBARFIT BRACE

FEATURES

- Anatomically-contoured arc stays fit lumbar curve and support lower back muscles.
- The narrow front design let user move freely without bulky feeling.
- Focused stays make brace great support on lumbar region.



SIZE	CM
S/M	65-95
L/XL	95-125

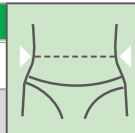

E-WA026


ADVANCED LUMBARFIT BRACE

FEATURES

- Anatomically-contoured arc stays fit lumbar curve and support lower back muscles.
- Additional front spiral springs restrict movement more effectively and provide superior support.
- Ventilated holes provide a comfortable fit and reduce heat and moisture.
- Extra side straps for fixing V-shaped belts add convenience of wearing.



SIZE	CM
S	65-80
M	75-90
L	85-100
XL	95-110
XXL	105-120

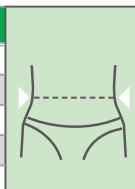

E-WA052

ABDOMINAL BINDER/

3 PANELS

FEATURES

- Provide compression to stabilize abdominal area
- Made of stretch elastic material
- Ideal for protection after surgery and support the belly



E-WA104

LUMBAR SACRO SUPPORT

FEATURES

- 4 rigid steel rods and 2 plastic rods support back muscles
- Effective tension side-pulls provide superior support
- Made of vented mesh material


E-WA106

LUMBAR SACRO SUPPORT

FEATURES

- 4 rigid anatomically contoured aluminum rods support back muscles
- Rhomboid tension side-pulls flatten the lordosis
- Breathable open cell foam material



SIZE	CM	
S	56-71	
M	69-84	
L	81-96	
XL	94-109	
XXL	107-122	

Suitable for E-WA052, E-WA104, E-WA106

E-WA113

DORSAL LUMBAR BACK SUPPORT

FEATURES

- Rigid aluminum uprights and pelvic bar limit spine flexion, extension, and lateral bending
- Bilateral shoulder straps with pads
- Made of foam padding material for maximal comfort



SIZE	CM	
S	56-71	
M	69-84	
L	81-96	
XL	94-109	
XXL	107-122	

Suitable for E-WA113, E-WA110

E-WA110

LUMBAR SACRO SUPPORT

FEATURES

- 4 rigid anatomically contoured aluminum rods support back muscles

- Effective vented tension side-pulls provide superior support
- Made of high quality vented elastic strapping material



E-WA120

OSTEOPOROSIS BRACE


FEATURES

- Anatomical aluminum plate fits spinal curve
- The front open piece makes this brace be easily applied like a vest
- Elastic bands constrain abnormal movement and posture but remain mobility
- Ideal for osteoporosis users



Aluminum upright



SIZE	CM	
S/M	56-84	
L/XL	84-109	

E-HE001E

HYPEREXTENSION BRACE


FEATURES

- Rigid frame restricts side-bending movement
- Moveable pelvic bar for better comfort while sitting
- Adjustable frame width and height
- Ideal for stable compression fractures, and osteoporosis



Quickly use buckle

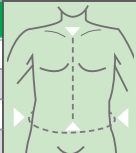

SIZE	PELVIC (CM)	LENGTH (CM)	
S	60-75	39-46	
M	75-90	42-49	
L	90-105	45-52	
XL	105-115	49-56	

E-HP008


HSO HIP SUPPORT

FEATURES

- 2 flexible and 3 rigid, removable aluminum stays help to control hip abduction/adduction
- Side pockets for convenient change of stays
- Wrap-around simple design
- Using of highly permeable fabric and skin-friendly material

SIZE	WAIST (CM)	THIGH (CM)
S	61-81	51-61
M	76-97	56-66
L	91-107	61-71
XL	102-127	66-81

LEFT /RIGHT

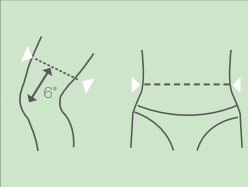


E-HP009

LITE HIP BRACE

FEATURES

- Rigid outer shell provides optimal support
- 6 blocking bolts limit hip flexion and extension
 - Extension limitation at: 0°, 20°, 40°, 60°, 80°, 100°
 - Flexion limitation at: 0°, 20°, 40°, 60°, 80°, 100°
- Abduction/ Adduction can be adjusted continuously between -60° to 60°
- Length adjustments from 42-48 cm
- Using of highly permeable fabric and skin-friendly material



SIZE
UNIVERSAL


E-HB002

HERNIA BELT

FEATURES

- Crotch straps are tied under the lower margin of gluteus maximum muscle to reduce excess pressure around buttock
- Two removable flexible compression cushions
- Made of velvet fabric material

SIZE	CM
S	63-76
M	76-89
L	89-102
XL	102-115

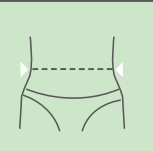




E-MB002

DELUXE MATERNITY BELT

FEATURES

- Supports the lumbar and abdominal regions to help relief of low back pain
- 2 rigid steel rods support back muscles
- Effective tension side-pulls provide superior support
- Made from lightweight tension elastic and velvet material

SIZE	WAIST CIRCUM.(CM)
S	Up to 103
M	103-113
L	113-123
XL	123-135

E-RB001

RIB BELT FOR MEN


E-RB002

RIB BELT FOR WOMEN

FEATURES

- Provide support to the rib cage
- Pain relief through compression whenever breathing or thorax movements are painful
- Made of high quality vented elastic strapping and velvet material
- Female breast contour design (E-RB002)

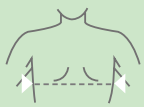

SIZE	E-RB002	E-RB001	 CHEST CIRCUM. (CM)
S	56-66	68-81	
M	66-76	81-96	
L	76-86	96-112	
XL	86-102	112-127	

NECK

E-CO007 / E-CO007B

ARMOUR CERVICAL COLLAR

FEATURES

- Two-piece polyethylene design enhances patient comfort and is easily adjusted with buckle design
- The unique adjustment technology allows 7 levels of chin height adjustment for customized fitting with a simple knob sliding.
- Removable heat pressing pads are made of soft for greater comfort and proper hygiene. The large tracheal opening provides functional access for the airway management and other essential procedures
- Safe for use with X-ray, CT, and MRI



PATENT NO.
CN206867320U
TWI605786B



1. Removable heat pressing pads ensure proper patient hygiene and comfort while preventing skin irritation



2. Buckle design makes applying easier and saves more time. Adjustment is only on the first applying by velcro.



3. The unique adjustment technology allows 7 levels of chin height adjustment for customized fitting with a simple knob sliding

E-CO007B

Velcro design makes hook and loop closure.



SIZE	NECK CIRCUM.	HEIGHT	
UNIVERSAL	30-45 CM	11-17.5 CM	


CO010 OB COLLAR

FEATURES

- OB collar is designed to provide great motion restriction without producing painful pressure points.
- Front and back open window design allow moisture and heat to be taken away.
- Trachea opening for advanced airway procedures, ventilation, drainage and sputum suction.
- Pad can be cleaned and reused.
- X-ray, CT and MRI compatible.



SIZE	CM
Pediatric	For 2~5 years old
S	<2
M	2-3
L	3-5
XL	>5





E-CO001/004 FOAM CERVICAL COLLAR/ FIRM DENSITY COLLAR

FEATURES

- Kinesthetic reminder to limit neck movement
- Made of foam with a stockinette cover
- Made of high density foam with a stockinette cover (E-CO004)

SIZE	CM
XS	25-32
S	32-38
M	38-43
L	43-48
XL	45-50

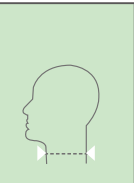

E-CO003


RIGID ADJUSTABLE COLLAR

FEATURES

- Anatomical shape and edge padding for great comfort
- Ventilation holes increase breathability
- Made of PE material with rigid support



SIZE	NECK (CM)	HEIGHT (CM)
S	30-35	7-8
M	35-40	8-9
L	40-45	9-10
XL	45-50	10-10.5

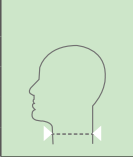

E-CO002A


CERVICAL COLLAR WITH TRACHEA OPENING

FEATURES

- Anatomical shape and preformed chin provide patient better comfort
- An anterior hole for a tracheostomy is available
- Made of closed cell foam material



SIZE	NECK (CM)	HEIGHT (CM)
S	25-33	5.7/8.3/ 10.8/13.3
M	33-41	
L	41-48	
XL	48 up	

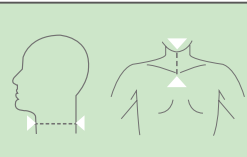

E-CH004

CERVICAL THORACIC ORTHOSIS

FEATURES

- Immobilize cervical spine and upper thoracic spine
- "Pull pin design" enables easy release for height adjustments of chin holder.
- All straps are adjustable
- Forehead strap adds stabilization
- Front and back supports for chin and occiput provide maximum control and suitable comfort



SIZE
UNIVERSAL

SHOULDER




E-AR002

SHOULDER IMMOBILIZER

FEATURES

- Padded cotton shoulder strap with velour-lined
- Envelope style with thumb support
- Breathable, durable material for arm rest

SIZE	CM
S	25-28
M	28-30
L	30-33
XL	33-36

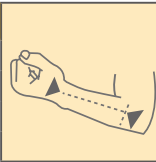



E-AR003

MESH ARM SLING

FEATURES

- Envelope style with thumb support
- Vented nylon mesh with shoulder strap

SIZE	CM
S	25-28
M	28-30
L	30-33
XL	33-36




E-AR007

SHOULDER SLING

FEATURES

- Wide chest swathe for optimal support
- Cross wearing provides better stabilization
- Comfortable, washable, skin-friendly terry toweling cloth material

SIZE
UNIVERSAL






E-AR011

SHOULDER IMMOBILIZER WITH ABDUCTION

FEATURES

- Immobilization of shoulder joint at 15° abduction
- Comfortable foam pillow provides more support
- An exercise ball facilitates hand function and stimulates circulation
- Quick release buckle makes it easy don and off
- Soft velour material

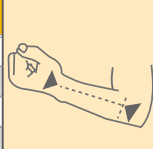SIZE	CM	
S	25-28	
M	28-30	
L	30-33	
XL	33-36	




E-CL003

CLAVICLE POSTURE SHOULDER BRACE

FEATURES

- Help to support an early kyphosis and to improve posture
- Flexible stays help provide added back support
- Made of high quality vented elastic strapping material

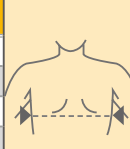SIZE	CM	
S	61-76	
M	76-91	
L	91-107	
XL	107-122	




E-CL004

CLAVICLE BRACE

FEATURES

- Help to alleviate clavicle fractures
- Excellent aid in correcting poor posture and shoulder slump
- Made of stretch elastic material for flexible compression

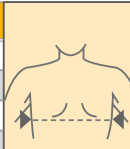
SIZE	CM	
S	61-76	
M	76-91	
L	91-107	
XL	107-122	

E-CL007

ADJUSTABLE POSTURE CORRECTOR

FEATURES

- The auxiliary straps on both sides can adjust the height of the shoulders
- Add two metal spring strips on the back to increase upper back support
- There are sponge pads in both categories to reduce the pressure and friction between the skin and the straps
- Breathable mesh material, more comfortable to wear

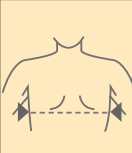SIZE	Height (CM)	Chest circumference (CM)	
S	Up to 155	Up to 65	
M	150-165	60-85	
L	165-180	75-100	
XL	175-190	90-115	




E-SH006

OVER SHOULDER HUMERAL SPLINT

FEATURES

- The two pieces of the brace come together to form a cylinder to provide compression.
- Shell of the deltoid brace is a rigid polyethylene that is lightweight and durable
- Brace is laminated with closed cell foam lining
- Deltoid cap extension of the broken shoulder brace prevents it from migrating
- The strap on proximal splint wraps beneath the opposite arm, limiting the motion of shoulder
- Ideal for fractured humeral treatment

SIZE	CM	
S	25-28	
M	28-30	
L	30-33	
XL	33-36	



Inner sleeve




E-SH001

NEOPRENE SHOULDER SUPPORT

FEATURES

- Elevate humeral head for anti-gravity
- Relief of tensile force at the subacromial space
- 1/8" thick neoprene with nylon lining
- Ideal for shoulder dislocation/subluxation

SIZE	CM	
S	25-28	
M	28-30	
L	30-33	
XL	33-36	
LEFT/RIGHT		



ELBOW

E-EL094A

POST-OP ELBOW BRACE WITH HAND GRIP

FEATURES

- Padded lined malleable aluminum cuffs for a perfect fit
- Simply press and rotate hinge with lock
- Extension limitation from 0° ~ 105°
- Flexion limitation from 0° ~ 120°
- Lock limitation from 0° ~ 90°
- Adjustable in every 15 degrees
- Anti-supination/pronation Arm Bar Kit is available



E-EL094

POST-OP ELBOW BRACE

FEATURES

- Provide single brace (without hand grip and strap) for customer's order



Padded lined malleable aluminum cuffs for a perfect fit

"Pull pin design" enables easy release for slide-to-size length adjustments.

Simply press and rotate hinge with lock

SIZE	LENGTH (CM)	CIRCUM. (CM)	
UNIVERSAL	34-43	Forearm Max. 32	
LEFT /RIGHT			

Suitable for E-EL094, E-EL094A



E-EL005

PNEUMATIC ELBOW BRACE

FEATURES

- The pneumatic pad provides relief at the muscle insertion point
- Hook and loop closure for pressure adjustment

SIZE	
UNIVERSAL	




E-EL006

PADDED TENNIS ELBOW BRACE

FEATURES

- The SBR pad provides compression at the muscle insertion point to relieve pain
- 1/8" thick neoprene with nylon lining

SIZE	CM	
S	23-25.5	
M	25.5-28	
L	28-30.5	
XL	30.5-33	




E-EL008

EPI

FEATURES

- Rigid plastic shell gives excellent support and provides compression
- Nonwoven fabric pad provides comfort

SIZE	CM	
XS	20-22	
S	22-25	
M	25-28	
L	28-31	
XL	31-34	

WRIST


E-WR123 FUNCTIONAL RESTING SPLINT

FEATURES

- Provide functional hand position, immobilization and prevent deformity.
- Optional finger separators.
- Prefabricated splint
- Plastic shell covered with foam padding cushion



SIZE	CM
S	6-7
M	7-8
L	8-9
XL	9-10
LEFT /RIGHT	




E-WR002 WRIST BRACE

FEATURES

- Provide compressive and proprioceptive effect over wrist joint.
- Wrap around design for precise fit
- Universal design to fit right and left hand


SIZE
UNIVERSAL

E-WR007 THUMB BRACE

FEATURES

- Immobilization of the thumb
- Anatomically-contoured aluminum thumb splint
- 1/8" thick neoprene with nylon lining

SIZE	CM
S	13-15
M	15-18
L	18-20
XL	20-22




E-WR050

WRISTLET WITH ABDUCTED THUMB SLEEVE

FEATURES

- Immobilization of the thumb
- Anatomically-contoured aluminum thumb splint
- Stretch elastic material for excellent compression
- Fit right and left hand



E-WR053

WRISTLET WITH FOREARM SPLINT

FEATURES

- Made of perforated suede with flannel lining for great comfort
- Removable malleable palmar stay for greater immobilization
- Stockinet liner provides comfort and easy applying
- Stretch elastic material for excellent compression



E-WR056

WRIST BRACE

FEATURES

- Immobilization and stabilization of the wrist
- Anatomically-contoured aluminum palmar splint
- Stretch elastic material for excellent compression

E-WR057

WRIST SPLINT

FEATURES

- Immobilization of the wrist
- Anatomically-contoured aluminum palmar splint
- V-notch at fifth MP for palmar spread
- Breathable open cell foam material





E-WR058
THUMB WRIST AND PALM SPLINT

FEATURES

- Immobilization of the wrist and thumb
- Anatomically-contoured aluminum palmar splint and thumb stay
- V-notch at fifth MP for palmar spread
- Breathable open cell foam material



E-WR060
WRIST BRACE WITH THUMB

FEATURES

- Anatomically-contoured aluminum thumb splint
- Extra-wide circumferential strap for superior support
- Breathable material with terry-lined interior

E-WR063
WRIST BRACE

FEATURES

- Anatomically-contoured aluminum palmar splint
- Extra-wide circumferential strap for superior support
- Universal design fits right or left hand
- Breathable material with terry-lined interior



SIZE	CM	
S	13-15	
M	15-18	
L	18-20	
XL	20-22	
LEFT /RIGHT		

Suitable for :
 WR050/ WR053/ WR056/ WR057/ WR058/ WR060/ WR063 (For both side)


E-WR071

SPICA THUMB WRIST LACER

FEATURES

- Anatomically contoured aluminum palmar splint
- Breathable – open cell structure, skin friendly material for comfort
- V-notch at fifth MP for palmar spread
- Easy single handed closure

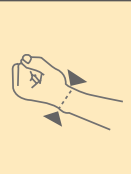

SIZE	CM	
S	14-16	
M	16-18	
L	18-20	
XL	20-22	
LEFT /RIGHT		


E-WR072

WRIST LACER

FEATURES

- Universal design for left and right hand
- Removable aluminum stays and anatomically contoured stay enhance stability
- Single pull lace closure system allows easy, one-handed application



SIZE	CM	
S/M	13-18	
L/XL	18-23	

E-FS001 BASEBALL FINGER SPLINT

FEATURES

- Immobilize sprained finger
- Stabilize the phalange joints in position
- Made of malleable aluminum with foam padding
- Packed in 3 pcs/set

SIZE	CM
S	8.9
M	10
L	12.7



E-FS004 FOUR PRONG FINGER SPLINT

FEATURES

- Immobilize sprained finger
- Stabilize the phalange joints in position
- Made of malleable aluminum with foam padding
- Packed in 3 pcs/set

SIZE	CM
S	3.8
M	5.0
L	7.6




E-FS002 TOAD FINGER SPLINT

FEATURES

- Immobilize sprained finger
- Stabilize the phalange joints in position
- Made of malleable aluminum with foam padding
- Packed in 3 pcs/set

SIZE	CM
S	5.4
M	7.0
L	9.5

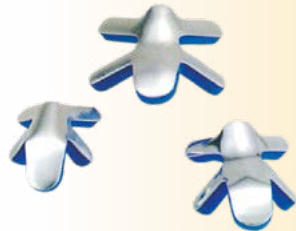

E-FS005 FOLD-OVER FINGER SPLINT


FEATURES

- Immobilize sprained finger
- Made of malleable aluminum with foam padding
- Packed in 3 pcs/set

SIZE	CM
S	6.4
M	7.0
L	8.9



SIZE	CM
S	6.4
M	10
L	15.2

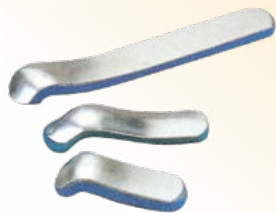

E-FS003 GUTTER FINGER SPLINT

FEATURES

- Immobilize sprained finger
- Stabilize the phalange joints in position
- Made of malleable aluminum with foam padding
- Packed in 3 pcs/set

E-FS006 MALLET FINGER SPLINT

FEATURES

- Stabilize the phalange joints in position
- Made of high quality rigid plastic.
- Ventilated design for great comfort.
- Packed in 3 pcs/set

SIZE
1
2
3
4
5
5.5
6
7





FOOTCARE

E-HV001

BUNION NIGHT SPLINT

FEATURES

- Plastic brace holds your toe in position to stretch out tight tendons and toe muscles.
- Protect bunion against rubbing
- Maintain metatarsal-phalangeal joint alignment



SIZE	SHOES (EU)
S	34-36
M	37-40
L	41-46
LEFT /RIGHT	

E-SI005

HEEL CUSHION

FEATURES

- Specially designed teardrop plug for additional heel cushion



SIZE	SHOES (EU)
S	34-38
M	38-42
L	42-46

Suitable for E-SI005, E-SI007

E-SI002

SILICONE HEEL CUP FOR SPURS-CENTRAL

FEATURES

- Cups shape provides heel support
- Specially designed soft blue dot for additional heel cushion



SIZE	SHOES (EU)
S	34-36
M	37-39
L	40-42
XL	43-45

E-SI007

"U" SHAPE HEEL PAD

FEATURES

- U-shaped design lifts up the heel and release the pain on the central area
- Made of medical level silicone rubber



E-SI021

FULL LENGTH SILICONE INSOLE

FEATURES

- Softer inserts absorb shock and pain on the metatarsal and heel areas.
- Raised metatarsal pad helps to support metatarsal arch.

SIZE	XS	S	M	L	XL	XXL	XXXL
SHOES (EU)	34	36	38	40	42	44	46



E-SI102

HEEL CUP WITH FABRIC

FEATURES

- Skin-friendly and hygienic fabric helps reduce friction
- Specially designed soft blue dot for local pressure and pain relief functions

SIZE	S	M	L	XL
SHOES (EU)	34-36	37-39	40-42	43-45



E-SI121

FULL LENGTH INSOLE WITH FABRIC

FEATURES

- Skin-friendly and hygienic fabric helps reduce friction
- Specially designed soft blue dot for local pressure and pain relief functions

SIZE	XS	S	M	L	XL	XXL	XXXL
SHOES (EU)	34	36	38	40	42	44	46





E-TS003

TOE SEPARATOR

FEATURES

- Made of medical level silicone rubber
- Reduce friction between toes

SIZE	S	M	L
------	---	---	---



E-TS005

BUNION SHIELD WITH LOOP

FEATURES

- Protect bunion against rubbing
- Maintain metatarsal-phalangeal joint alignment

SIZE	S/M	L/XL
------	-----	------



E-TS009

TUBULAR TOE BANDAGE

FEATURES

- Absorb shock and shear force to reduce pain
- Made of TPR

SIZE	LENGTH
XS	15 CM
S	
M	
L	



E-TS014

FOOT –BUNION SEPARATOR

FEATURES

- Absorbs lateral shock and minimizes friction on the bunion area.
- It brings constant relief of the affected site.
- Recommend S/M size for 5th toe bunion and L/XL for big toe.

SIZE	S/M	L/XL
------	-----	------

PEDIATRIC





E-EL881

ELBOW IMMOBILIZER

FEATURES

- To immobilize or stabilize your child's elbow after an injury, sprain, bone fracture, or surgery.
- Made from a comfortable, latex free foam material that is soft against the skin.

SIZE	SPLINT HEIGHT (CM)	ELBOW CIRCUM.(CM)	
PEDIATRIC	10	Up to 19	
YOUTH	15	18-23	



E-WR801

WRIST SPLINT

FEATURES

- Anatomically-contoured aluminum palmar splint
- Low contour at palmar crease allows full finger function
- Made of comfortable, washable, skin-friendly terry toweling cloth material

SIZE	AGE IN YEARS	WRIST CIRCUM.(CM)	
PEDIATRIC	2-7	11-14	
YOUTH	8-12	14-17	
LEFT /RIGHT			



E-WR802

WRIST SPLINT WITH ABDUCTED THUMB

FEATURES

- Rigid thumb stay for custom fit
- Anatomically-contoured aluminum palmar splint
- Made of comfortable, washable, skin-friendly terry toweling cloth material

SIZE	AGE IN YEARS	WRIST CIRCUM.(CM)	
PEDIATRIC	2-7	11-14	
YOUTH	8-12	14-17	
LEFT /RIGHT			

E-KN851

NEOPRENE HINGED KNEE SUPPORT

FEATURES

- Open patellar design
- Hinges hidden in pockets offer superior stability
- Made of closed cell neoprene




E-KN870

NEOPRENE STABILIZED KNEE SUPPORT

FEATURES

- Open patellar design
- Hinges hidden in pockets offer superior stability
- Made of closed cell neoprene



SIZE	AGE IN YEARS	KNEE CIRCUM.(CM)	
PEDIATRIC	6 - 7	25 – 30	
YOUTH	8-12	30 – 33	

Suitable for : E-KN851, E-KN870

E-KN881

KNEE IMMOBILIZER

FEATURES

- Removable medial, lateral and posterior stays immobilize the leg
- Available length: 8", 12"
- Made of comfortable, washable, skin-friendly terry toweling cloth material with pediatric print



E-KN882

KNEE IMMOBILIZER/ 3 PANELS

FEATURES

- Three piece design for precise fit
- Removable medial, lateral and posterior stays
- Available length: 8", 12"
- Made of soft velvet material with washable and skin-friendly terry-lined interior



SIZE	LENGTH	KNEE CIRCUM. (CM)	
UNIVERSAL	8", 12"	25 – 30	

Suitable for : E-KN881, E-KN882




E-CO801

SOFT FOAM COLLAR

FEATURES

- Kinesthetic reminder to limit neck movement
- Contour design for optimal fit
- Made of medium density foam with a stockinette cover

SIZE	NECK CIRCUM.(CM)	HEIGHT(CM)	
YOUTH	25-31	6	




E-CO802A

CERVICAL COLLAR WITH TRACHEA OPENING

FEATURES

- Restrict cervical spine flexion & extension with gentle support
- Light weight, with trachea opening
- Kinesthetic reminder to limit neck movement
- Made of closed cell foam material, prevent shear stress

SIZE	NECK CIRCUM.(CM)	HEIGHT(CM)	
PEDIATRIC	25-33	6	



E-CO806

ADJUSTABLE COLLAR

FEATURES

- Anterior oversized trachea opening for emergency
- Back openings for ongoing cervical assessment
- Radiolucent, MRI and CT scan compatible
- Made of polypropylene material with smooth foam padding

SIZE	NECK CIRCUM.(CM)	
PEDIATRIC	20-46	

E-HM801 CRANIAL PROTECTION HELMET

FEATURES

- Made of material IXPE with excellent impact resistant
- Article ventilation isolation provides great ventilation between head and hat
- Y type double loops reduce the migration of the hat during activities
- Adjustable head circumference is able to fit most
- Two colors available (blue and pink)




SIZE	HEAD CIRCUM.	
UNIVERSAL	43cm~55cm	

E-AR801 ARM SLING WITH PAD

FEATURES

- Soft shoulder strap with cotton pads
- Envelope style with thumb support
- Made of comfortable, washable, skin-friendly terry fabric

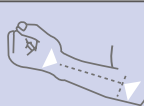


SIZE	AGE IN YEARS	FOREARM LENGTH (CM)	
PEDIATRIC	2-7	13-20	
YOUTH	8-12	18-25	



E-CL801 CLAVICLE SUPPORT

FEATURES

- Help to correct and remind poor posture and Immobilize fractured clavicle
- Using of foam strap padding with stockinet covered

SIZE	AGE IN YEARS	CHEST CIRCUM.	
PEDIATRIC	2-7	51-61	

E-CL803

ADJUSTABLE CLAVICLE SUPPORT

FEATURES

- Help to correct and remind poor posture and Immobilize fractured clavicle
- Size adjustable design
- Made of stretch elastic material for flexible compression

SIZE Universal



E-CL804

POSTURE CORRECTOR

FEATURES

- Aluminum alloy stays provide optimal support
- Adjustable support level by axilla straps
- Made of breathable material
- Comfortable shoulder pads



Two strap fixed ring provides different support level.



SIZE	CM	
S	30-40	
M	40-50	
L	50-65	
XL	65-80	

E-WK802B

WALKING BOOT

FEATURES

- Contoured strut design provides better support and fit
- Rocker bottom design assists normal gait
- Washable foam liner improves fit and comfort
- Front and back panels for superior protection (E-WK802)



E-WK802



SIZE	CHILD SHOE(US)
S	<5-6
M	7-8
L	9-10
XL	11-12

E-AN860

AIR / GEL STIRRUP ANKLE BRACE

FEATURES

- Anatomically-contoured plastic shells control inversion or eversion
- Combined air and gel bladders for cold therapy and even-distributed pressure
- Foam padding with tricot liner provides good shock-absorption



SIZE	AGE IN YEARS	ANKLE CIRCUM.	
PEDIATRIC	2-7	13-19	




E-CA801

CAST BOOT

FEATURES

- Rocker bottom design assists normal gait
- Non-skid sole
- Durable canvas material with adjustable closures

SIZE	LENGTH(CM)	WIDTH(CM)	
PEDIATRIC	18	9	
YOUTH	22	10	




E-PH801

PAVLIK HARNESS

FEATURES

- Harness is a dynamic flexion abduction orthosis
- Booties design provides advanced comfort
- Soft stirrups with contact closure for ease of application and adjustment

SIZE	MONTHS	CHEST(CM)	
XS	0-3	30-36	
S		36-41	
M	3-6	41-46	
L	6-9	46-53	
XL	9-12	53-60	




E-PH802

FREJKA SPLINT

FEATURES

- Ideal for congenital hip luxation and sub-luxation, hip dysplasia.
- Flexible polyethylene splint with foam-padded edges for comfort covered by nice textile diaper

SIZE	(CM)	
XS	14	
S	16	
M	18	
L	20	
XL	22	
XXL	24	

Tied on both sides of strap

eLife®
embrace life.





eLife

SPORTS

E-KND502

POWERGUARD KNEE PAD

FEATURES

- High impact 3.5cm thickness foam pad for optimal impact absorption
- Elastic band with anti-slip designed for preventing slipping down
- Stretch lycra mesh for long lasting comfort

SIZE	CM
UNIVERSAL	38-52



Measure 10cm above knee joint



FEATURES

- Carbon fiber frame with contoured shape and rigid frame for excellent fitting and stability
- Flexible subshell with adjustable strap for more customized fit
- Polycentric hinges with plug-in flex/ext stops
Extension limitation at: 0°, 10°, 20°, 30°, 40°
Flexion limitation at: 0°, 10°, 20°, 30°, 45°, 60°, 75°, 90°
Immobilization limitation at: 0°, 10°, 20°, 30°
- Numbered straps provide desirable counter force and adjust optimal size.
- Breathable padding liner for superior ventilation and skin contact

SIZE	THIGH(CM)	CALF(CM)
S	39-47	29-37
M	47-54	37-43
L	54-60	43-50
XL	60-67	50-57

LEFT /RIGHT



E-KNK035

X1 KNEE BRACE



E-KNK035A

X1 KNEE BRACE WITH CUP





E-ELD502

**POWERGUARD
ELBOW BRACE**

FEATURES

- High impact EVA pad for superb impact absorption
- Elastic band with anti-slip designed for preventing slipping down
- Stretch lycra mesh for long lasting comfort

SIZE	CM
S	20-25
M	24-29
L	28-33



E-KND501

**POWERGUARD
KNEE BRACE**

FEATURES

- High impact EVA pad for superb impact absorption
- Elastic band with anti-slip designed for preventing slipping down
- Stretch lycra mesh for long lasting comfort

SIZE	CM
S	30-37
M	36-43
L	42-50



PROTECTION



PROTECTION >>>

PROTECTION

elfe



E-ELS270
PRESSINGBAND ELBOW
BRACE

FEATURES

- Banded compression improve muscular endurance and protect the elbow and relieve muscle soreness, pain and reduce swelling
- Circular knitting technology with high quality materials to produce the elbow area compression

SIZE	CM
S	19-22
M	22-25
L	25-28
XL	28-31



E-SNS270
PRESSINGBAND CALF
BRACE

FEATURES

- Banded compression improve muscular endurance and protect the calf and relieve muscle soreness, pain and reduce swelling
- Circular knitting technology with high quality materials to produce the calf area compression


SIZE	CM
S	30-34
M	34-38
L	38-42
XL	42-46





PRESSINGBAND

TPU partial pressurized to improve muscular endurance and protect the extremities, reduces inflammation and swelling. Annulus patterns design provides additional support and compression on peripheral joints and muscle tissues.

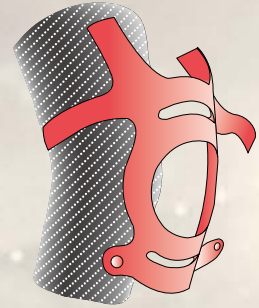

E-ANS270 **PRESSINGBAND** **ANKLE BRACE**

FEATURES

- Banded compression improve muscular endurance and protect the ankle and relieve muscle soreness, pain and reduce swelling
- Circular knitting technology with high quality materials to produce the ankle area compression

SIZE	CM
S	16-21
M	21-26
L	26-31
XL	31-37



E-KNS270 **PRESSINGBAND** **KNEE BRACE**

FEATURES

- Banded compression improve muscular endurance and protect the knee and relieve muscle soreness, pain and reduce swelling
- Circular knitting technology with high quality materials to produce the knee area compression

SIZE	CM
S	38-43
M	41-45
L	44-49
XL	47-52



KNITTING

E-life Sports Medicine Knitting Brace are developed for general sport activities. PressingBand series provide constant and stable joints and soft tissues support that improve exercise performance. E-life also have Flat knitting product. Flat knitting technology makes braces provide optimal body fit. E-life Sports Medicine Knitting Brace is ideal for sporty support while engaging in active mobility.

E-KNS453

CAMOUFLAGE KNEE BRACE

FEATURES

- Three-dimensional flat knitting technology, achieve the partial conversion from normal muscle pressure
- Flexible spiral stays in both sides avoid the support to be out of shape
- Patella silicone buttress can increase the stability of the knee

SIZE	CM
S	33-36
M	36-38
L	38-41
XL	41-46



E-KNS252

REINFORCED KNIT KNEE BRACE

FEATURES

- Reinforced straps provide additional support
- Flexible spiral stays in both sides avoid the support to be out of shape
- Patella silicone buttress can increase the stability of the knee

SIZE	CM
S	38-43
M	41-45
L	44-49
XL	47-52





E-ANCO01
COOL-FIT
ANKLE
SUPPORT

FEATURES

- Ideal for mild to moderate ankle distortions and ankle instability
- Cross type straps provide extra compression
- Hook and loop closure design fits most ankle sizes

SIZE

UNIVERSAL



E-ELC001
COOL-FIT
ELBOW
SUPPORT

FEATURES

- Ideal for chronic inflammation, general protection
- Provides comfortable compression and reduces pain
- Hook and loop closure design fits most elbow sizes

SIZE

UNIVERSAL



E-WRC001
COOL-FIT
WRIST
SUPPORT

FEATURES

- Ideal for slight wrist strain or sprain
- Wrap-around design for super fit

SIZE

UNIVERSAL





SIZE
UNIVERSAL



E-KNC001
COOL-FIT
WRAP KNEE SUPPORT

FEATURES

- Hook and loop closure design fits most knee sizes
- Low pressure of patella design
- Ideal for mild muscle strain around knee joint and articular swelling



SIZE	CM
S/M	30-38
L/XL	38-46



E-KNC051
COOL-FIT
KNEE BRACE

FEATURES

- Front open Velcro type is simple fitting
- Bilateral contour spiral stays provide better stability
- Ideal for reducing stress on tendons and muscles around knee

COOL-FIT



E-KNC058
COOL-FIT
WRAP KNEE SUPPORT

FEATURES

- Front open Velcro type is simple fitting
- Bilateral contour stays provide better stability
- Ideal for reducing stress on tendons and muscles around knee

SIZE	CM
S/M	39-43
L/XL	42-46



E-KNC072
COOL-FIT
WRAP AROUND KNEE BRACE

FEATURES

- Wrap around design provides easy application
- Hook and loop closures above and below knee for secure fit
- Ideal for reducing stress on tendons and muscles around knee

Air mesh



SIZE	CM
S/M	39-43
L/XL	42-46



COOL-FIT >>>

COOL-FIT BRACE



elife



E-ELN002
THERMO ELBOW BRACE

- Increase blood circulation through warm effect
- Contoured design for freedom of movement

SIZE	CM
S	22-25
M	25-28
L	28-31
XL	31-34



E-ELN003
THERMO ELBOW SUPPORT

- Increase blood circulation through warm effect
- Contoured design for freedom of movement
- Additional adjustable tension straps

SIZE	CM
S	22-25
M	25-28
L	28-31
XL	31-34



E-SHN002
THERMO SHOULDER BRACE

- Increase blood circulation through warm effect
- Contoured design for freedom of movement

SIZE	CM
S	80-90
M	90-100
L	100-110



E-SNN002
THERMO CALF BRACE

- Increase blood circulation through warm effect
- Contoured design for freedom of movement

SIZE	CM
S	30-34
M	34-38
L	38-42
XL	42-46



E-THN002
THERMO THIGH BRACE

- Increase blood circulation through warm effect
- Contoured design for freedom of movement

SIZE	CM
S	44-49
M	49-54
L	54-59
XL	59-65





E-WRN002
THERMO WRIST BRACE

- Increase blood circulation through warm effect
- Contoured design for freedom of movement
- Additional adjustable tension straps

SIZE	CM
S	13-16
M	16-19
L	19-22
XL	22-25



E-WAN002
THERMO LUMBAR BRACE

- Increase blood circulation through warm effect
- Contoured design for freedom of movement

SIZE	CM
S	70-80
M	80-90
L	90-100
XL	100-110
2XL	110+



E-WAN003
THERMO LUMBAR SUPPORT

- Increase blood circulation through warm effect
- Contoured design for freedom of movement
- Flexible stays provide additional lumbar support
- Four plastic stays provide additional lumbar support

SIZE	CM
S	7080
M	80-90
L	90-100
XL	100-110



E-ANN002
THERMO ANKLE BRACE

- Increase blood circulation through warm effect
- Contoured design for freedom of movement

SIZE	CM
S	19-22
M	22-25
L	25-28
XL	28-31



E-ANN003
THERMO ANKLE SUPPORT

- Increase blood circulation through warm effect
- Contoured design for freedom of movement
- Flexible stays for enhanced support

SIZE	CM
S	19-22
M	22-25
L	25-28





E-KNN002
THERMO KNEE BRACE

FEATURES

- Increase blood circulation through warm effect
- Contoured design for freedom of movement
- Open patellar design

SIZE	CM
S	30-35
M	35-40
L	40-45
XL	45-50



E-KNN003
THERMO KNEE SUPPORT

FEATURES

- Increase blood circulation through warm effect
- Contoured design for freedom of movement
- Open patellar design
- Flexible stays in both sides avoid the support to be out of shape

SIZE	CM
S	30-35
M	35-40
L	40-45
XL	45-50



E-KNN051
THERMO HINGED KNEE SUPPORT

FEATURES

- Increase blood circulation through warm effect
- Contoured design for freedom of movement
- Open patellar design
- Pivot hinge in both sides avoid the support to be out of shape

SIZE	CM
S	30-35
M	35-40
L	40-45
XL	45-50




E-EPN005
THERMO SLIM SHORTS TONER

FEATURES

- Increase blood circulation through warm effect
- Contoured design for freedom of movement

SIZE	CM
S	70-80
M	80-90
L	90-100
XL	100-110



THERMO SUPPORT 

**THERMO
SUPPORT**



E-life Sports Medicine is committed to developing the ideal brace that offers complete protection for athletes and sport lovers. Braces varying in strength are available to offer all-round protection of injured joints, muscles, ligaments, tendons, and other tissues, thereby ensuring excellent stability.

To cater to the needs in sports arenas or outdoor activity environments, E-life Sports Medicine made of neoprene material produces warmth and moderate pressure on

joints and muscles. In addition, cool-fit material has also been developed to provide good ventilation and moisture absorption as well as environmental protection material features. Knitting technology is also adopted to provide moderate pressure and the 3D trim for a more comfortable fit.

Enthusiastically adhering to the R&D spirit and innovative ideas, we strive to create more benefits for sport doers in the field of sports medicine.

SPORTS MEDICINE

THERMO SUPPORT	103
COOL-FIT BRACE	99
KNITTING	96
PROTECTION	93



elite



eLife®

2020

ORTHOPEDICS

enrich
life.